Colamark















Colamark A33 High Speed Labeling & Plunger Assembly System for Prefilled Syringes

Designed for use in the packaging line of pre-filled syringes, the system labels on the prefilled syringe and at the same time assemble the plunger-rod into the syringe. The DN1Denester can be used for automatic syringe feeding from from the transfer nest. Optional inspection station and assembly station for plunger backstopper.



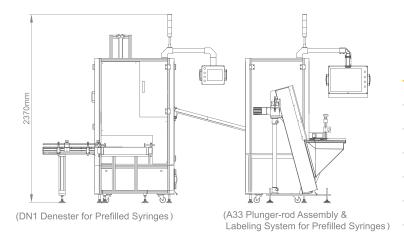

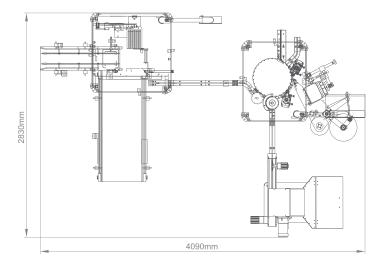

Model

Features

- **Syringes** unloaded from full nest. plunger-rod fed by vibrating feeder, indexing turntable with stations for inserting plunger-rod, tightening, and labeling.
- Robust design and use of world-class materials and components to ensure exceptional performance, stability and quality.
- The labeling system uses the proven emery traction roll and servo drive technology to ensure extra precision and long term stability.
- Various on-line coding options, including hotstamping, inkjet, laser, thermal transfer are available to meet customer's specific needs.
- Various inspection options are available to detect missing label, missing code, wrong code, missing plunger-rod etc. to ensure quality output for the user.
- Defective items are automatically collected in a rejection tray.

DN1 Denester for Prefilled Syringes A33 Plunger-rod Assembly & Labeling System for Prefilled Syringes

Model	, , , , , , , , , , , , , , , , , , , ,
Material size	Diameter: 8.15±0.1mm / 10.85±0.1mm
System capacity	Up to 300 ppm (depending on material size)
Labeling precision	±0.3-0.5mm (excluding errors from materials)
Denesting method And speed	±0.3-0.5mm (Robotic operated denesting comb. Denested syringes are fed into the in-feeding chute of A33. Capacity up to 300 ppm.errors from materials)
Compressed air	0.4-0.7Mpa
Electricity	220VAC±10% 46-62Hz 2.0kW
Dimension	(L)4,090 × (W)2,830 × (H)2,370 mm

(Custom engineered solutions are available to cater for customer's special needs.)



150 William Street Perth Amboy, NJ 08861 1-800-298-8936 Phone: 732-324-0800 Fax: 732-324-4488 Email: sales@dartronics.com www.dartronics.com

Colamark (Guangzhou) Labeling Equipment Limited

19 Jin Xiu Road, GETDD, Guangzhou, China

Tel: +86 20 8221 6718 Fax: +86 20 8221 5799 PC:510730

Colamark Technologies Limited 820 Topsail Plaza,11 On Sum Street, Shatin, Hong Kong Tel: +852 2634 8011 Fax: +852 8161 3311 Email: info@colamark.com Web: www.colamark.com

For continuous improvement on design and product innovations, all specification data are subject to change without prior notice

| HongKong | Guangzhou | Kunming Chengdu | Wuhan | Chongqing | Xian | Shanghai | Qingdao | Beijing | Shenyang | Harbin |